

INSTRUCTION MANUAL

Printer Terminal





This is a hazard alert mark.



This mark informs you about the operation of the product.

Note This manual is subject to change without notice at any time to improve the product. No part of this manual may be photocopied, reproduced, or translated into another language without the prior written consent of the A&D Company.

Product specifications are subject to change without any obligation on the part of the manufacture.

Copyright@1999 A&D Company, Limited

Correction Sheet

This is correction sheet for the printer terminal instruction manual of TM-2480 and TM-2481, version 1600-1A-IE.

Page	Incorrect	Correct
r aye		
42	10.3.4 Alarm Setting for Measurement Time	10.3.4 Alarm Setting for Measurement Time
	Item number Setting value	Item number Setting value
42	Item number	tem number 1-1 1-2 2-1 2-1 2-2 3-1 3-2 4-1 4-2 5-1 5-2 6-1 6-2
43	Setting value Press D 1 10 SSEED D 2 0 SSEED D 3 3 0	Item number Setting value I - I II I - Z II

Page	Incorrect		Correct
9		Symbols	
M. M. Carlos		ψ	Turning on / off the recorder
			Direct current.
No.		~	Alternating current.
		1999	Date of manufacture.
		Δ	Attention symbol. "See instruction for use".
			Class II equipment.
		C €0366	CE marking.
n.b.		SN	Serial number.

Note:


The following function is disabled on the TM-2480/TM-2481 Printer Terminal.

On the UA PC series (MODEL 2), the alarm setting function is disabled.

Related pages:

Chap. 10.3 Model 2, Operation and Setting - page 37 (shown below)

Chap. 10.3.4 Alarm Settings for Measurement Time - page 42

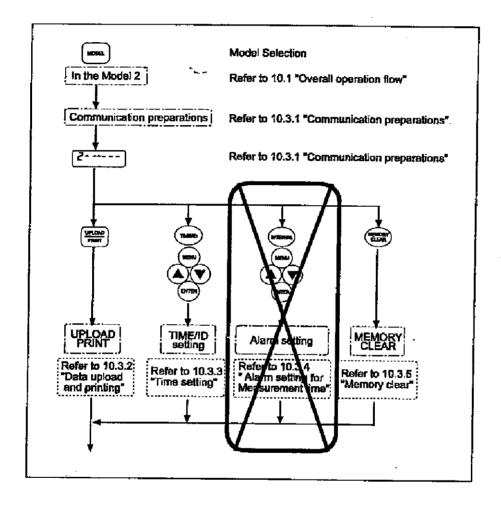

Contents

1.	Compliance 3
1.1.	Compliance with European Directive 93/42 EEC for Medical Product · · · 3
1.2.	Compliance with FCC Rules 3
2.	Safety 3
3.	Packing list····· 4
4.	Introduction · · · · · 5
5.	Features · · · · · 5
6.	Specifications 7
7.	Part Names · · · · · 8
8.	Explanation of Key Switches 10
9.	Preparation · · · · · · · · · · · · · · · · · · ·
9.1.	Cable Connection 11
9.2.	Power (AC adapter or batteries)
9.3.	Printer Paper Installation (Replacement) · · · · · 12
10.	The Whole Operation · · · · · 13
10.1.	The Whole Operation
10.1.1	. Model Selection · · · · · · 16
10.2.	Model 1, Operation and Setting · · · · · 18
10.2.1	. Communication Preparations 20
10.2.2	2. Data Upload and Printing · · · · · 22
10.2.3	B. Print Format Setting 24
10.2.4	Time/ID Setting 26
10.2.5	6. Measurement Interval Setting · · · · · 28
10.2	2.5.1. Setting of Automatic Measurement Mode III
10.2.6	Memory Clear 35
10.2.7	Printing of Setting Content
10.3.	Model 2, Operation and Setting · · · · · · 36
10.3.1	. Communication Preparations 37
10.3.2	Data Upload and Printing

10.3.3.	Time Setting · · · · · · 40
10.3.4.	Alarm Setting for Measurement Time 42
10.3.5,	Memory Clear · · · · · · 44
10.3.6.	Printing of Setting List · · · · 45
10.4.	Model 3, Operation and Setting · · · · · 46
10.4.1.	Communication Preparations · · · · · 47
10.4.2.	Data Upload and Printing · · · · · 48
11. M	aintenance · · · · · · 49
11.1,	Cleaning 49
11.2.	Before Requesting Repairs · · · · · 49
	ptions····· 50
13. A	ppendix ····· 50
13.1.	External Diagrams · · · · 50
13.2.	Printout Examples 52

1. Compliance

1.1. Compliance with European Directive 93/42 EEC for Medical Product

The device conforms to European Directive 93/42 EEC for Medical Products. This is evidenced by the CE mark of conformity accompanied by the reference number of a designated authority.

1.2. Compliance with FCC Rules

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. (FCC = Federal Communications Commission in the U.S.A.)

2. Safety



Printer Terminal

 Be sure to read the instruction manual of the connected blood pressure monitor before use.

Batteries

- Use alkaline batteries in the printer. Use of other types of batteries may cause reduced printer performance, damage, or fire.
- Do not mix and use old and new batteries.
- When the printer will not be used for a long period of time, remove the batteries.
 Batteries kept in the printer for a long period may leak and cause damage.

Storage Location

- Do not store the printer in:
- Locations where water or other fluids may get on or into the printer.
- Locations with high temperature or humidity.
- Locations in direct sunlight.
- Locations where the printer will be subject to excessive vibration or shocks.
- Locations with dust, salt, sulfur, or other contaminants.

Repair

Do not allow the case to be opened by anyone other than authorized service personnel. This will void the warranty and may cause damage to the printer or fires. When repairs are necessary, contact the retailer where you purchased this product or yours nearest A & D representative.

Cleaning

When the printer terminal requires cleaning, gently wipe with a cloth dampened with water or a neutral cleanser.

3. Packing list

Package Content

Confirm that the following items are in the package.

Printer terminal

Thermal paper: 1 roll
Communication cable: 1
AC adapter: 1
Instruction manual: 1

4. Introduction

Thank you for purchasing the Printer Terminal (TM-2480/TM-2481) from A & D. This manual was written for the TM-2480/TM-2481. Please read the manual carefully before using the unit in order to ensure a sufficient understanding of the printer terminal and provide proper use. Store the manual in a readily available location.

5. Features

- The printer can be used both ambulatory blood pressure monitors and blood pressure monitors used in the home.
- Blood pressure measurement results can easily be printed.
- The TM-2481 has an infrared communication function.
- When an ambulatory blood pressure monitor with infrared communication functions is used, measurement data can be transferred without a cable (TM-2481 only).
- The settings of clock and automatic measurement functions can also be input in blood pressure monitors.
- Applicable blood pressure monitors of the TM-2430 series.
 The TM-2430 series is "The portable automatic blood pressure monitors".
 - Every 24 hours the unit statistically processes blood pressure, pulse, and pressure load data and then prints the results.
 - Every 24 hours the blood pressure, pulse, and pressure load data are printed as trend graphs.
 - Statistical processing is done to calculate the average value and standard deviation for systolic pressure, diastolic pressure, and pulse rate. The results are then printed.
 - The pressure load is calculated for systolic and diastolic pressure using a surface area that exceeds a freely set value. (Units: mmHg•h)
 - The basic period for statistical processing is 24 hours. Results can also be printed for interval blocks and freely set waking and sleeping periods.

- Applicable blood pressure monitors of the TM-2431 series.
 The TM-2431 is "The portable blood pressure monitor with infrared communication function".
 - Data can be read and settings made through infrared communications without a cable connection (TM-2481 only).
- Applicable blood pressure monitors of the UA PC series and UB PC series.
 The UA PC series is "Principally upper arm blood pressure monitor for family use" and the UB PC series is "Wrist blood pressure monitor for family use".
 - Every month the unit statistically processes blood pressure data and pulse data and then prints the results.
 - Every month the blood pressure data and pulse data are printed as trend graphs.
 - Statistical processing is done to calculate the average value and standard deviation for systolic pressure, diastolic pressure, and pulse rate. The results are then printed.
 - The number of times the systolic/diastolic pressure exceeds 140/90 is counted. When measurement data for two months or more is present, all blood pressure data is printed. Then all blood pressure trend graphs are printed. Finally, all pulse data trend graphs are printed.

6. Specifications

Power

AC adapter or 6 x alkaline batteries

(type LR6, type AA, Mignon)

Power consumption

9.6W max.

Display

LCD

Printer

Serial thermal printer

Print speed

Approx. 1line/sec.

Printer paper

50mm diameter X 58mm width

Communication

functions

Serial communications:

EIA RS232C standard based Infrared communications: IrDA SIR1.0 standard based

(TM2481 only)

Operating

temperature and humidity ranges

+10 ~ +40°C 85%RH or less

(however, there must be no condensation)

Storage temperature

-10 ~ +55°C

and humidity ranges

95%RH (however, there must be no condensation)

External dimensions

TM-2480: Approx. 220mm (W) X 130mm (D) X 85mm (H)

TM-2481: Approx. 220mm (W) X 180mm (D) X 90mm (H)

Weight


TM-2480: 530g

TM-2481: 580g

Accessories


Printer paper, communication cable, AC adapter

7. Part Names

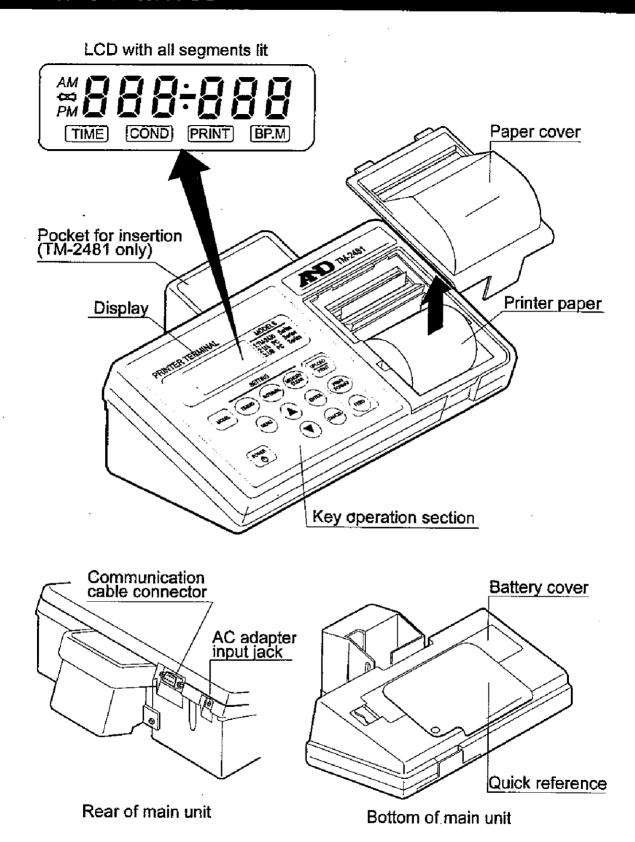

Connected Model Display

When the power switch is turned on, the LCD will blink and then enter the standard mode, as shown below. The number at the left of the LCD at this time shows the model number of the blood pressure monitor. The correlation between model numbers and model names is shown in the table below.

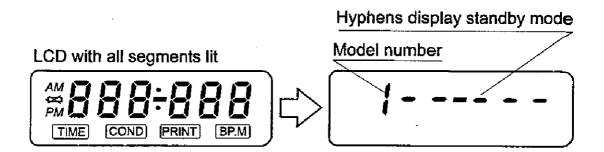


Correlation between Model Numbers and Model Names

Model number	Model name
1	TM-2430 Series
2	UA PC Series
3	UB PC Series

Explanation of Display Marks

LCD with all segments lit

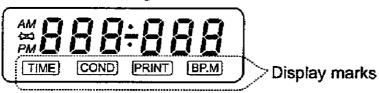


Display mark	Explanation
TIME	Displayed when setting the time or ID (applicable to model numbers 1 and 2).
COND	Displayed when setting the measurement interval (applicable to model numbers 1 and 2).
PRINT	Displayed when setting the report format (applicable to model number 1).
-PRINT)-	Blinks during printing.
BP.M	Displayed when measurement data has been uploaded.


Note: Unless this is a back light display, it can not flash.

Low Battery Indicator

When the battery voltage drops during use, the "Lb" is displayed. Replace new batteries or use the AC adapter.

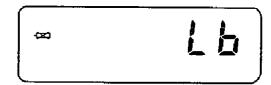

8. Explanation of Key Switches

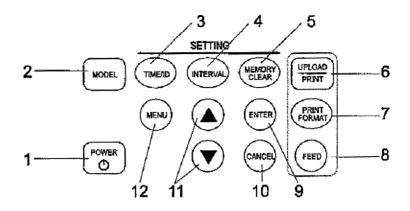

Number	Key	Explanation
1	POWER	Turns the power on and off.
2	MODEL	Selects the connected blood pressure monitor model.
3	TIMEND	Enters the mode for setting time and ID.
4	INTERVAL	Enters the mode for setting the measurement interval.
5	MEMORY	Cléars the measurement data memory of the blood pressure monitor.
6	UPLOAD PRINT	Uploads measurement data from the blood pressure monitor. Printing is then done using the report format in memory.
7	PRINT	Used to enter the mode for setting the report format.
8	FEED	Feeds printer paper.
9	ENTER	Inputs setting changes when in a setting mode.
10	CANCEL	Cancels communications, printing, or settings.
11		Changes setting values when in a setting mode.
12	MENU	Selects the item to be changed when in a setting mode.

9. Preparation

Perform the following preparations before using the printer.


9.1. Cable Connection

If you use the infrared communication, it is not necessary to connect the cable.

Step 1 Insert the communication cable into the communication cable connector on the back of the main printer unit, and secure the cable.

Notice

- When connecting the UB PC Series, the cable is different. In this case, use the special communications adapter to connect the cable.
- There are some shapes regarding the accessory cable.

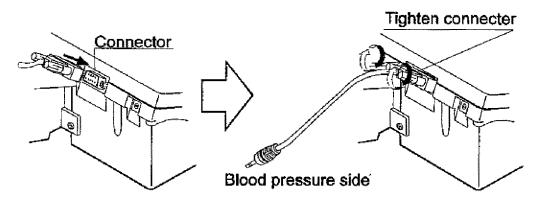


9.2. Power (AC adapter or batteries)

This printer can use batteries or the AC adapter.


Using the AC Adapter

When using the AC adapter, refer to the diagram at right and insert the accessory AC adapter into the AC adapter input jack on the back of the printer. Plug the other end of the adapter into a wall socket.

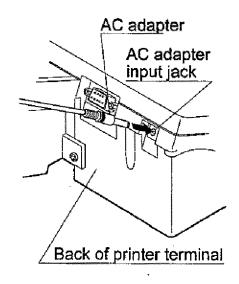

Using Batteries

When using batteries, turn the printer over, refer to the diagram at right, and remove the battery holder from the printer terminal. Follow the arrows in the diagram to put in the batteries.

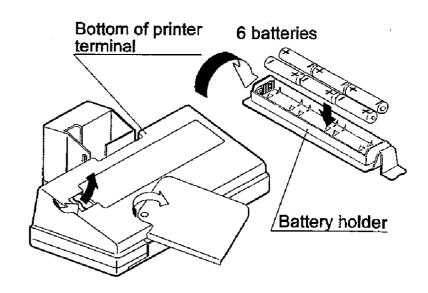

⚠ Caution!

- Put in the batteries properly according to the + and signs displayed on the holder.
- Replace all batteries when replacement is necessary.
- Do not use different types of batteries together.
- When the printer will not be used for a long time, remove the batteries to avoid damage from leakage.
- When using batteries, the printer terminal is automatically turned off by the auto-power-off function.

9.3. Printer Paper Installation (Replacement)


Notice

- This printer uses thermal paper. Use only the specified printer paper.
- When only 60cm of paper remains, the pink end marker will appear. Replace the paper at this time.


Printer Paper Installation Method

Step 1 Remove the paper cover.


Step 2 Cut the end of the printer paper as shown in the diagram at right.

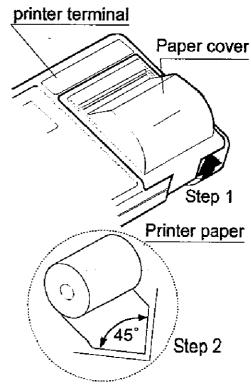

Step 3 Insert the end into the specified position.

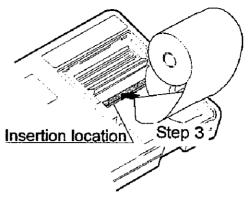

Step 4 When the end comes out of the printer, pull slowly.

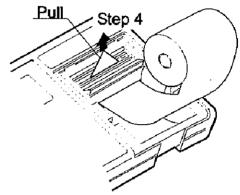

- Step 5 When the entire cut portion of the end is visible, press the FEED key and confirm that the paper is fed properly.

 If the paper is not fed properly, repeat the process from step 2.
- Step 6 If the paper feeds properly, replace the paper cover over the printer paper and secure.

10. The Whole Operation

10.1. The Whole Operation

The overall operation flow is shown here for model 1 (TM-2430 series). Be sure to properly understand this flow when using the units.

Although the flow is the same for other units, refer to the operation and setting sections for a detailed explanation.

- Step 1. Connect the cable


 Refer to section 9.1 " Cable Connection".
- Step 2. Connect the AC adapter. (Or install the batteries.)
 Refer to section 9.2 "Power" for details.
- Step 3. Install the printer paper Refer to section 9.3 " Printer Paper Installation" for details.
- Step 4. Press the MODEL key and select the model being connected.

Model number	Model name
1	TM-2430 series
2	UA PC series
3	UB PC series

Refer to section 10.1.1 "Model selection" for details.

- Step 5. Press the TIME/ID key and set the time.
- Step 6. Press the TIME/ID key to set the clock. Since the time and ID can be set at the same time, it is recommended that changes be made to both values simultaneously. Refer to section 10.2.4 "Time/ID setting" for details.
- Step 7. Press the INTERVAL key and set the automatic measurement interval for the blood pressure monitor.

 Refer to section 10.2.5 " Measurement Interval Setting " for details.

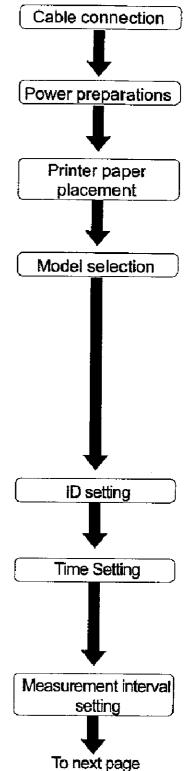

- Step 8. Press the ENTER key. Printing of the blood pressure monitor and printer terminal setting will begin.


 Refer to section 10.2.7 " Printing of Setting Content " for details.
- Step 9. Measure the patient's blood pressure.

 Refer to the operation manual of the blood pressure monitor for details.
- Step 10. Press the UPLOAD/PRINT key. The data will be uploaded and printing will automatically begin.

 Printing will be done using the previously set report format.
- Step 11. When the uploaded data is to be printed in a different report format, press the PRINT FORMAT key and select the desired format.

 Refer to section 10.2.3 " Print Format Setting " for details.
- Step 12. The data recorded in memory should be deleted in preparation for the next blood pressure measurement. Press and hold the MEMORY CLEAR key for 5 seconds or more to clear the data.

 Refer to section 10.2.6 "Memory Clear " for details.

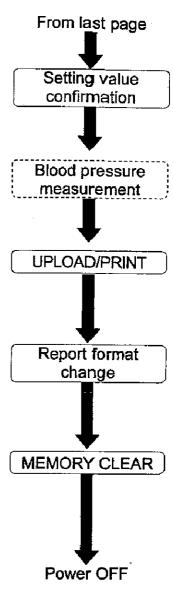


Notice

- Operation or setting that are not required can be skipped.
- When setting do not require changes, data processing is possible using only "Step 10 " and "Step 12 ".

10.1.1. Model Selection

When the overall operation flow is sufficiently understood and printer preparations are complete (items in section 9. "Preparation" are complete), first select the blood pressure monitor model to be connected.

- Step 1. Press the POWER key. When the power is turned on the unit will enter the standby mode.
- Step 2. Press the MODEL key and select the model of the connected blood pressure monitor. The relationship between model numbers and blood pressure monitor model names is shown below.

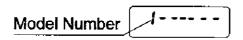


Model number	Model name	Reference section
1	TM-2430 series	10.2 Model 1 operation and setting content
2	UA PC series	10.3 Model 2 operation and setting content
3	UB PC series	10.4 Model 3 operation

Subsequent operation

Once the model number of the connected blood pressure monitor is selected, the following operation will be executed. However, there are some differences depending on the model. Use the reference sections listed in the above table as necessary.

- Step 3. Connect the blood pressure monitor to the printer (except when using infrared communications) to prepare for communication between the units.
- Step 4. Upload measurement data, set the various setting, and select the report format.
- Step 5. When each operation is complete, the unit will return to the standby mode. The operation in step 4 may be repeated at this point if necessary.
- Step 6. When operation is complete, press the POWER key to turn the power off.

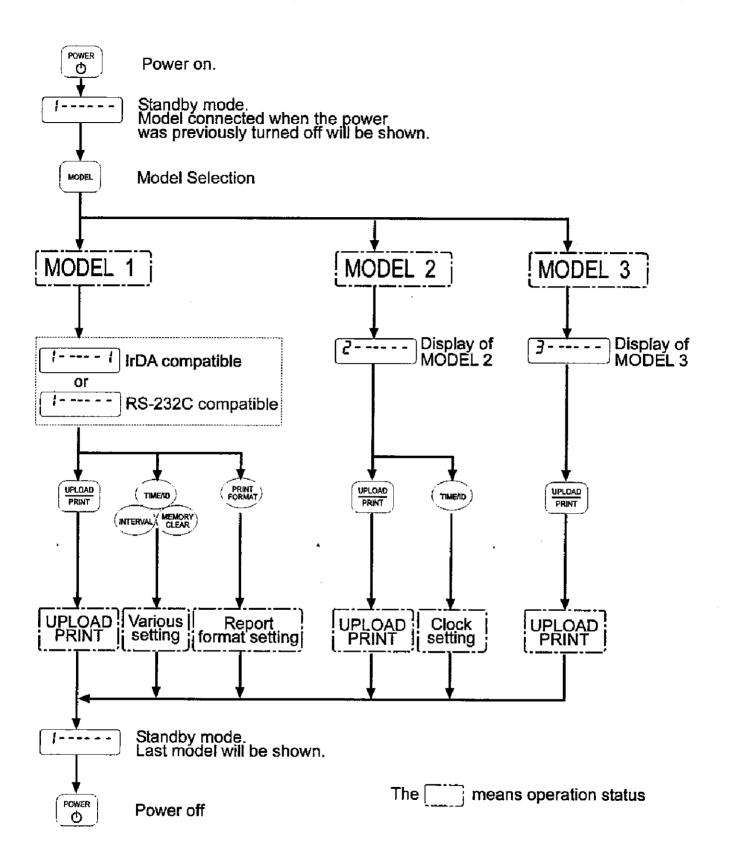

10.2. Model 1, Operation and Setting

When model 1 is selected (10.1.1 "Model selection"), operation is outlined in the following flowchart.

Model name: TM-2430 series Turning on -Standby mode. Model connected when the power was previously turned off will be shown. MODEL Model Selection In the Model 1 Refer to 10.1 "Overall operation flow" Communication preparations Refer to 10.2.1 "Communication preparations" Communication method selection IrDA compatible Refer to 10.2.1 "Communication preparations" or RS-232C compatible UPLOAD PRINT FORMAT TIMETO (INTERVAL) PRINT MENU MENU MENU ENTER ENTER) ENTER) Report UPLOAD TIME/ID INTERVAL MEMORY format setting setting <u>selection</u> Refer to 10.2.3: Refer to 10.2.2 Refer to 10.2.4 Refer to 10.2.5 Refer to 10.2.6 "Data upload "Print format "Time / ID "Measurement "Memory clear setting" and printing" setting" interval setting": In the Model 1 means operation status.

The various relevant setting used when model 1 is selected are shown below. With model 1, all of the printer terminal (TM-2480/2481) functions can be used. For a detailed description of the setting, refer to the reference section listed in the table below. In addition, with model 1 infrared (TM-2481) can be selected as the communication method. Refer to section 10.2.1 "Communication preparations" for details regarding these selections.

Function name	Content	Reference section
UPLOAD/PRINT	Used to upload and print measurement data	10.2.2 "Data upload and printing"
PRINT FORMAT setting	Sets the report format for the measurement data	10.2.3 "Print format setting"
TIME/ID setting	Sets the time and ID number. Display of the blood pressure monitor measurement result can also be turned on or off.	10.2.4 "Time / ID setting"
INTERVAL setting	Sets the interval for automatic measurement by the blood pressure monitor.	10.2.5 "Measurement interval setting"
MEMORY CLEAR	Clears the measurement data in the memory of the blood pressure monitor.	10.2.6 "Memory clear"

10.2.1. Communication Preparations


With model 1, cable communication or infrared communication (TM-2481 only) can be selected. A brief description of these communication methods is provided below.

Cable Communication (RS232C communication)

- Connect the accessory communication cable to the connector on the back of the Step 1. printer terminal (refer to section 9. "Preparation").
- Press the **A** ∇ keys to select the communication method. Step 2. When the far right digit on the LCD display is a hyphen (-), cable communication is selected.



- Connect the mini-pin jack end of the accessory communication cable to the Step 3. blood pressure monitor (TM-2430 series). The blood pressure monitor (TM-2430 series) display will change to [----]
- When the display changes, preparations for cable communication are complete. Step 4.

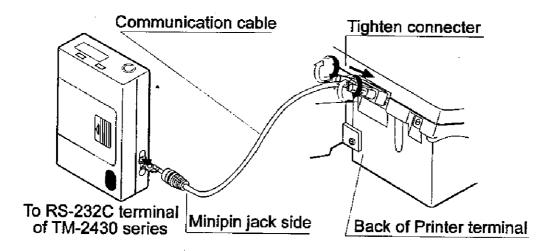

Caution!

When the unit is left with the communication cable inserted in the blood pressure monitor, the batteries will rapidly deteriorate. Therefore, when communications are complete, unplug the cable.

When the cable is unplugged the blood pressure monitor will return to the clock display.

Infrared Communications (IrDA communication) - TM-2481 only


Notice

Infrared communications are limited to use with blood pressure monitors having the infrared communication function (TM-2431).

Step 1. Switch the printer terminal to infrared communications mode.

Press the keys to select the communication method.


When the far right digit on the LCD display is a 1, infrared communication is selected.

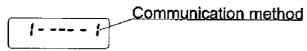

- Step 2. Insert the blood pressure monitor (TM-2431) into the back panel pocket of the TM-2481.
- Step 3. Changing the blood pressure monitor (TM-2431) to infrared communication mode

Press the AUTO ON/OFF switch on the blood pressure monitor (TM-2431).

- Step 4. If the switch is held down the blood pressure monitor display will soon read
- Step 5. When the display changes, preparations for infrared communications are complete.

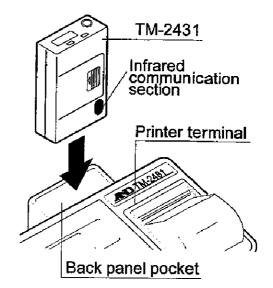

Λ

Caution!

When the blood pressure monitor is left in the infrared communications mode (in the infrared communications are complete, press the AUTO ON/OFF switch to exit the infrared communication mode.

When the infrared communications mode is exited the blood pressure monitor will return to the clock display.

10.2.2. Data Upload and Printing

Explanation

With model 1, the print format can be selected. The various print formats are described in the table below.

Print format	Content	
1	Prints a table + blood pressure graph + pulse graph + pressure load graph (refer to section 13. Appendix (TM-2430 series printout example)).	
2	Prints a table + blood pressure graph + pulse graph.	
3	Prints a table + pressure load graph.	
4	Prints a table only.	
5	Prints the blood pressure graph only.	
6	Prints the pulse graph only.	
7	Prints the pressure load graph only.	


- The method for setting the report format, including the print formats for model 1, is described in section 10.2.3 "Print format setting."
- When the UPLOAD/PRINT key is pressed, measurement data from the blood pressure monitor is uploaded. The measurement results are then printed according to the selected report format.
- The report format is saved even when the power is turned off.
- When measurement data is uploaded, the BPM indicator is displayed. While this mark is displayed the printer will maintain the measurement data in memory. If the report format is changed while the measurement data is recorded in memory, printing will begin using the new report format.
- Uploaded measurement data will be deleted when the connected model is changed or the power is turned off.

Setting Procedure


- Step 1. Confirm that communication preparations are complete (refer to section 10.2.1 "Communication preparations").
- Step 2. Press the UPLOAD/PRINT key.

 Measurement data will be uploaded from the blood pressure monitor.
- Step 3. The BPM indicator is displayed when uploading is complete.
- Step 4. Printing will begin. The PRINT mark will blink during printing.
- Step 5. The PRINT indicator will go out when printing is complete.


Example of cable communication

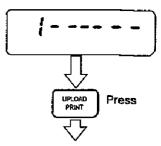


After measurement datauploading is complete

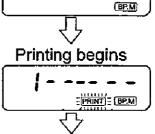

Display after printing is complete.

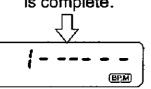

10.2.3. Print Format Setting

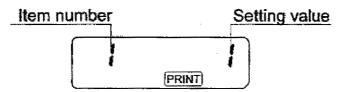

Explanation


- The report format setting is used to input various parameters into the printer terminal. This setting does not change the measurement data, time, or other setting of the blood pressure monitor.
- The above display appears when the PRINT FORMAT key is pressed.
- The PRINT indicator displayed during report format setting.
- When measurement data is in memory (BP.M) indicator is displayed), printing in the new format will begin after changes are made.
- The report format is saved in memory even when the power is turned off.
- An explanation of item numbers and setting values is provided in the table below.

Item	Setting	Catting contact	
number	value range	Setting content	
•	1~7	Print format 1 Table + blood pressure graph + pulse graph + pressure load graph 2 Table + blood pressure graph + pulse graph 3 Table + pressure load graph 4 Table 5 Blood pressure graph 6 Pulse graph	
2	AU, 0 ~ 23	7 Pressure load graph Processing start time AU Time of first measurement data (automatic) 0 ~ 23 Set time	
3	0 ~ 23	Sleep time	
4	0 ~ 23	Awake time	
5	100 ~ 180	Systolic threshold value while awake (set in 5mmHg steps)	
6	50 ~ 130	Diastolic threshold value while awake (set in 5mmHg steps)	
7	100 ~ 180	Systolic threshold value while a sleep (set in 5mmHg steps)	
8	50 ~ 130	Diastolic threshold value while a sleep (set in 5mmHg steps)	
9		Initial setting Print format: Processing start time: Sleep time: AU Sleep time: AWake time: To (7:00am) Systolic threshold value while awake: Diastolic threshold value while awake: Systolic threshold value while a sleep: Diastolic threshold value while a sleep:	


Setting Procedure

- Step 1. When communication preparations are complete (refer to section 10.2.1 "Communication preparations"), press the PRINT FORMAT key.
- Step 2. The PRINT indicator is displayed and item number 1 will be displayed along with the current setting value.
- Step 3. To change the setting value, press the ▲ ▼ keys.
 In this example the setting value in item number 1 (print format) will be changed to 3 (table + pressure load graph).
- Step 4. When changing the item number, press the MENU key. When changing the setting value, press the ▲ ▼ keys.

 In this example the setting value of item number 2 (processing start time) will be changed to 3 (3:00AM).
- Step 5. When changes are necessary, repeat the same basic operations to change the setting values for each item number.
- Step 6. When changes are complete press the ENTER key to confirm the new values. When measurement data is (BP.M) contained memory in indicator displayed), printing will begin using the new setting. When measurement data is not contained in memory, "PLEASE UPLOAD DATA" will be printed and the report format will be set.
- Step 7. If the CANCEL key is pressed without pressing the ENTER key, report format setting will be canceled and the report format will not be changed.

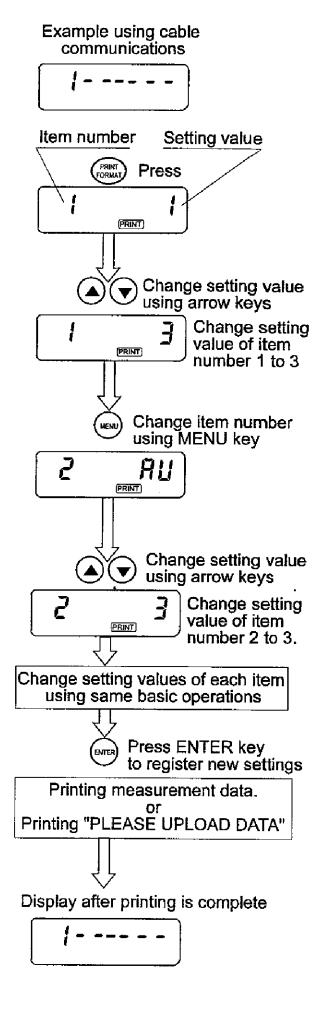

10.2.4. Time/ID Setting

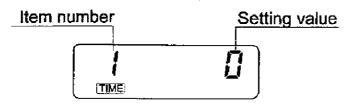

Explanation

- The time/ID setting sets the various parameters in the blood pressure monitor (TM-2430 series).
- The display shown above appears when the TIME/ID key is pressed.
- The TIME indicator is displayed during time/ID setting operation.
- An explanation of item numbers and setting values is shown in the table below.

Item number	Setting value range	Setting content
1	0,1	Display function on/off 0: Off , 1: On
2	1 ~ 99	ID number
5	1997 ~ 2096	Year (4 digits)
6	1 ~ 12	Month
7	1 ~ 31	Day (setting range is automatically adjusted according to the year and month)
8.	0 ~ 23	Hour (use 24-hour time)
9	00 ~ 59	Minutes (2 digits are always displayed)


Use the display function on/off item only after reading the operation manual of the blood pressure monitor.

Setting Procedure


- Step 1. When communication preparations are complete (refer to section 10.2.1 "Communication preparations"), press the TIME/ID key.
- Step 2. The TIME indicator and item number 1 will be displayed along with the setting value for the connected blood pressure monitor.
- Step 3. To change the setting value, press the ▲ ▼ keys.
 In this example the setting value in item number 1 (display function on/off) will be changed to 1 (on).
- Step 4. When changing the item number, press the MENU key. When changing the setting value, press the ▲ ▼ keys.

 In this example the setting value of item number 2 (ID number) will be changed to 90.
- Step 5. When other changes are necessary, select each item number and repeat the same basic operations to change the setting values for year, month, day, hour, and minute.
- Step 6. When changes are complete press the ENTER key to confirm the new values.

 The new setting for the blood pressure monitor will be printed. Please check the content of the setting.
- Step 7. If the CANCEL key is pressed without pressing the ENTER key, time/ID setting operation will be canceled and the setting of the blood pressure monitor will not be changed.

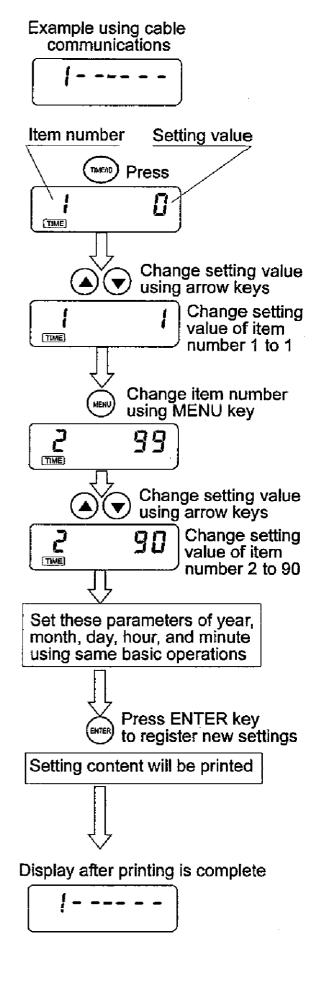

10.2.5. Measurement Interval Setting

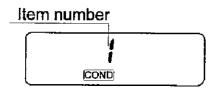

Explanation

- The measurement interval setting sets the measurement interval to be used for automatic measurement by the blood pressure monitor (TM-2430 series).
- The above display appears when the INTERVAL key is pressed.
- The COND indicator is displayed during measurement interval setting.
- An explanation of the item numbers is provided in the table below.

Item number	Setting content	
1	Automatic measurement mode I 07:00 ~ 21:59: 15-minute intervals 22:00 ~ 6:59: 30-minute intervals	
2	Automatic measurement mode II 15-minute intervals when S is not displayed 30-minute intervals when S is displayed	
3	Automatic measurement mode III Maximum of 6 blocks can be set Enter other setting routine	


For details regarding the measurement interval refer to the operation manual of the blood pressure monitor (TM-2430 series).

Setting Procedure

- Step 1. When communication preparations are complete (refer to section 10.2.1 "Communication preparations"), press the INTERVAL key.
- Step 2. The COND indicator is displayed and the item number set in the connected blood pressure monitor will be displayed.
- Step 3. To change the item number, press the MENU key.
 In this example the item number 1 (automatic measurement mode I) will be changed to 2 (automatic measurement mode II).
- Step 4. When changes are complete press, the ENTER key to confirm the new value.


 (When item number 3 is selected, setting operation will continue. Refer to the next section for automatic measurement mode III setting.)

 The new setting for the blood pressure monitor will be printed. Please check the content of the setting.
- Step 5. If the CANCEL key is pressed without pressing the ENTER key, the measurement interval setting operation will be canceled and the setting of the blood pressure monitor will not be changed.

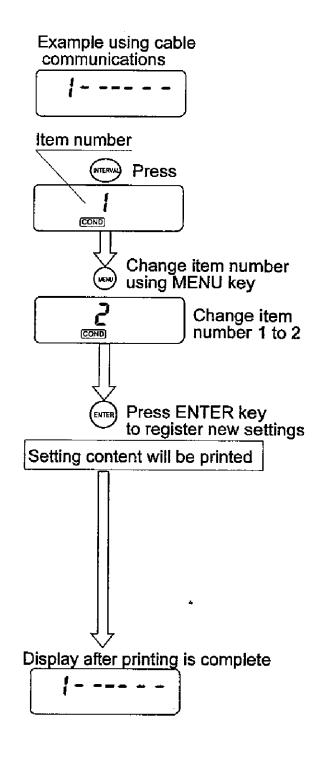

10.2.5.1. Setting of Automatic Measurement Mode III

This is automatic measurement mode III (item number 3) setting.

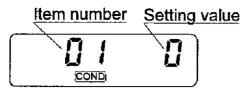

Explanation

- When automatic measurement mode III is selected as the measurement interval, the 24-hour time period in a day can be divided into up to 6 blocks and the measurement interval can be set for each block.
- The COND indicator will continue to be displayed and the above display will appear.
- An explanation of each item number is provided in the table below.

Item number	Setting range	Setting content
01	0 ~ 23	Block 1 start time
02	(Note 1)	Measurement interval
03	0 ~ 23	Block 1 end time (block 2 start time).
04	(Note 1)	Measurement interval
05	0 ~ 23	Block 2 end time (block 3 start time). However, a value within the previous block time cannot be set.
06	(Note 1)	Measurement interval
07	0 ~ 23	Block 3 end time (block 4 start time). However, a value within the previous block time cannot be set.
08	(Note 1)	Measurement interval
09	0 ~ 23	Block 4 end time (block 5 start time). However, a value within the previous block time cannot be set.
10	(Note 1)	Measurement interval
11	0 ~ 23	Block 5 end time (block 6 start time). However, a value within the previous block time cannot be set.
12	(Note 1)	Measurement interval
13	Block 1 start time	Block 6 end time However, the block 1 start time must be input.

Note 1: The selective range is OFF, 5, 10, 15, 20, 30, 60, and 120 in units of minute.


Setting Procedure

- Step 1. Enter the mode 1 using the MODEL key form the standby mode. (Refer to section10.2.1).
- Select Step 2. item the 3 as interval measurement and press the ENTER key. Setting automatic measurement mode III will then be possible. The setting method for block 1 time. start measurement interval, and end time will be used as an example for this explanation. The setting method for subsequent blocks is the same.
- Step 3. Item number 01 (block 1 start time) will be displayed. When changing the start time setting value, press the ▲ ▼ keys. In this example the setting value of item number 01 (block 1 start time) is changed to 1 (start time of 1:00AM).
- item number 02 (block 1 measurement interval). When changing the measurement interval setting value, press the ▲ ▼ keys.

 In this example the setting value of item number 02 (block 1 measurement interval) is changed to 30 (30 minutes).

Press the MENU key to display

Step 5. Press the MENU key to display item number 03 (block 1 end time). When changing the end time setting, press the A V keys. The block 1 end time automatically becomes the block 2 start time. In this example the setting value of item number 03 (block 1 end time) is changed to 22 (10:00PM).

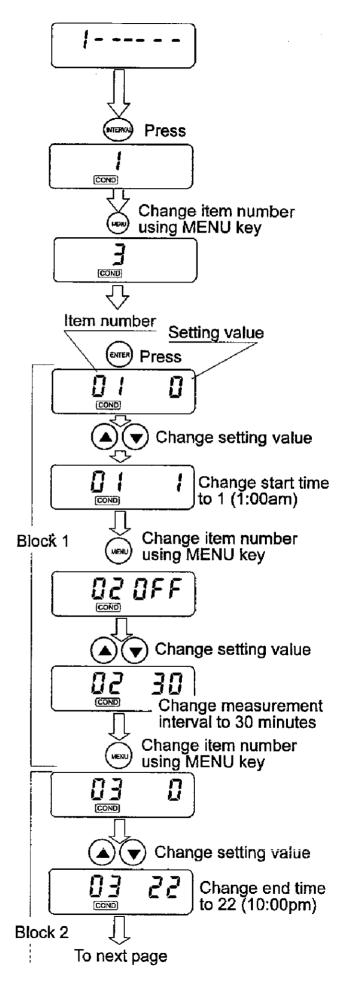


Step 4.

Step 6. When setting for block 2 and subsequent blocks necessary. use the same procedure select to item numbers and change setting values. However, be sure to carefully read and understand the setting cautions described in the next section.

Step 7. When changes are complete, press the ENTER key to register the new values and exit the setting mode.

The new setting for the blood pressure monitor will be printed. Please check the content of these setting.

Step 8. If the CANCEL key is pressed without pressing the ENTER key, the measurement interval setting operation will be canceled and the setting of the blood pressure monitor will not be changed.

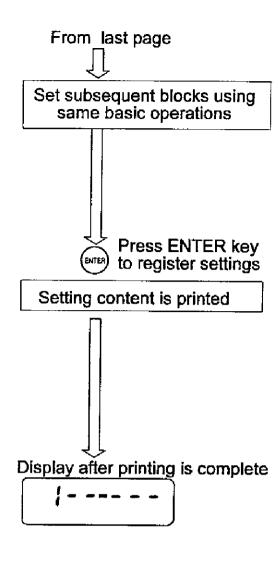


Setting Cautions

• The time range set in a previous block cannot be set in the next block.

Example

Block 1: 1 ~ 5 (1:00AM ~ 5:00AM) Block 2: 5 ~ 14 (5:00PM ~ 2:00PM)

In this case the setting value for block 3 end time (item number 07) must be within the $14 \sim 23$, $0 \sim 1$ range.

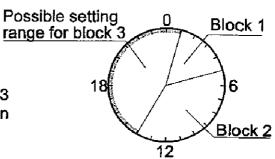

Item	Setting	Setting content	
number	range		
01	1	Block 1 start time	
02		Measurement interval	
03	5	Block 1 end time (block 2 start time)	
04	-	Measurement interval	
05	14	Block 2 end time (block 3 start time)	
06	_	Measurement interval	
07	14 ~ 23 0 ~ 1	Block 3 end time (block 4 start time). However, the time range used in previous blocks cannot be set.	

 Even when a total of 6 blocks are not used, further block setting is not possible once a total of 24 hours has been used by all currently set blocks.


Example

Block 1: 1 ~ 5 (1:00AM ~ 5:00AM)

Block 2: 5 ~ 14 (5:00AM ~ 2:00PM)

Block 3: 14 ~ 1 (2:00PM ~ 1:00AM)

In this case the block 4 measurement interval setting (item number 08) will not appear when the menu key is pressed. Instead, item number 01 will appear.

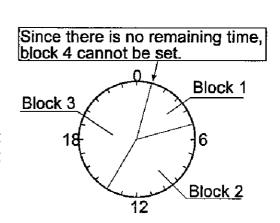

 When 120 minutes is selected as the measurement interval, the display will require an end time using 2-hour intervals.

Example

Block 1 start time: 1 (1:00AM)

Block 1 measurement interval: 120 minutes

In this case the setting values displayed for block 1 end time (item number 03) will be 1, 3, 5, 7, 9, 11, 13, 15, 17, 19, 21, and 23.

 The 120 minute measurement interval will not be displayed for final blocks containing only 1 hour of time.

A final block is defined as the next block when the total time set for previous blocks is 23 hours.


Example

Block 1 start time: 0 (12:00AM)

Block 1 measurement interval (60 minutes)

Block 1 end time 23 (11:00PM)

In this case, a setting value of 120 cannot be input for the measurement interval of block 2 (item number 04).

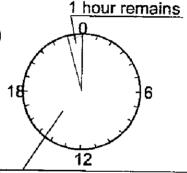

Block 1

The situation is the same when multiple blocks use this period of time.

10.2.6. Memory Clear

Explanation

The MEMORY CLEAR key is used to delete the measurement data in the blood pressure monitor (TM-2430 series).

Setting Procedure


Step 1. When communication preparations are complete (refer to section 10.2.1 "Communication preparations"), press and hold the MEMORY CLEAR key. The buzzer will beep every second and, after 5 seconds, a long buzzer sound will be heard.

After measurement data has been cleared, the message "MEMORY CLEARED" will be printed.

When measurement data has not been cleared, the message "MEMORY CLEAR FAILED" will be printed.

Step 2. When the MEMORY CLEAR key is not pressed and held for 5 seconds the measurement data will not be deleted.

In this case no message will be printed

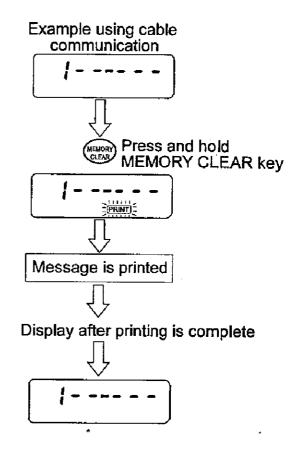

⚠ Caution!

Be careful when using these operations since measurement data cleared from the blood pressure monitor cannot be recovered. Be sure to execute measurement data processing before clearing the memory.

When attaching the blood pressure monitor to a new person, be sure to clear the measurement data.

10.2.7. Printing of Setting Content

Explanation


- The time, measurement-interval, ID, and other items set in the blood pressure monitor (TM-2430 series) will be printed.
- The report format set in the printer terminal (TM-2480) will also be printed.

Setting Procedure

Step 1. When communication preparations are complete (refer to section 10.2.1 "Communication preparations"), press the ENTER key. Printing will begin.

Notice


When no blood pressure monitor is connected, only the report format will be printed.

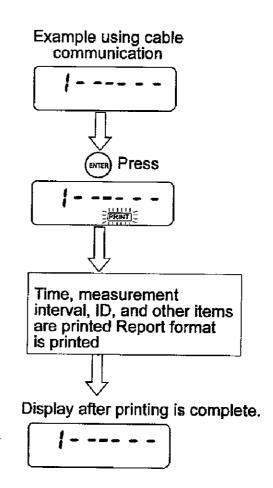

10.3. Model 2, Operation and Setting

When model 2 is selected in the overall operation flow (10.1.1 Model selection) shown in section 10.1, the subsequent operations are outlined in the following flowchart.

Model name: UA PC series

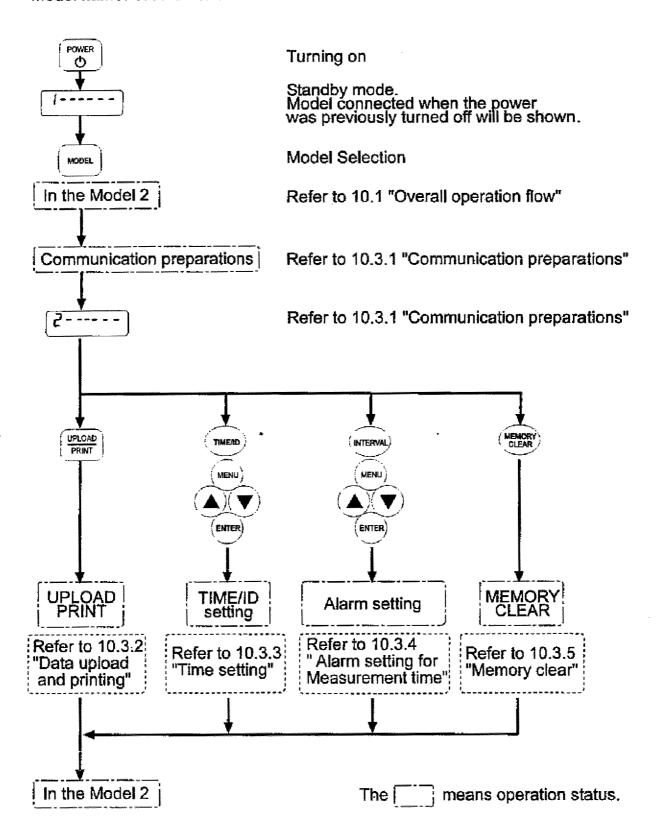


When model 2 is selected the functions that can be used with the printer terminal (TM-2480) are limited. The functions and various setting that can be used are shown in the table below. For a detailed description of the setting, refer to the reference section listed in the table below.

Function name	Content	Reference section
UPLOAD/PRINT	Used to upload and print measurement data	10.3.2 "Data Upload and Printing"
Time Setting	Sets the time	10.3.3 "Time Setting"
Alarm Setting	Sets the alarm for measurement of the blood pressure monitor	10.3.4 " Alarm Setting for Measurement Time"
Memory Clear	Clears the measurement data in the blood pressure monitor	10.3.5 "Memory Clear"

10.3.1. Communication Preparations

- Step 1. Connect the accessory communication cable to the connector on the back of the printer terminal (refer to section 9. "Preparation").
- Step 2. When the blood pressure monitor displays clock, connect the mini-pin jack of the accessory communication cable to the blood pressure monitor.

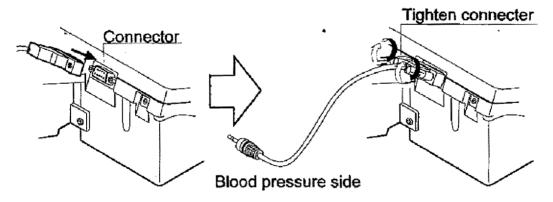

⚠ Caution!

When the unit is left with the communication cable inserted in the blood pressure monitor, the batteries will rapidly deteriorate. Therefore, when communications are complete, remove the cable. When the cable is removed, the blood pressure monitor will return to the clock display.

About several products, this connection procedure is different. Refer to the instruction manual of product you use

10.3.2. Data Upload and Printing

Explanation


- When the UPLOAD PRINT key is pressed, measurement data is uploaded from the blood pressure monitor. Measurement results are then printed using the standard report format.
 - With model 2, the report format cannot be selected.
- When measurement data is uploaded, the BP.M indicator is displayed.

 Measurement data uploading is complete when this indicator is displayed. The communication cable can be unplugged from the blood pressure monitor at this time.


Setting Procedure

- Step 1. Confirm that communication preparations are complete (refer to section 10.3.1 "Communication preparations").
- Step 2. Press the UPLOAD/PRINT key.

 Measurement data will be uploaded from the blood pressure monitor.
- Step 3. The BP.M indicator is displayed when uploading is complete.
- Step 4. Printing will begin. The PRINT indicator will blink during printing.
- Step 5. The PRINT indicator will go out when printing is complete.

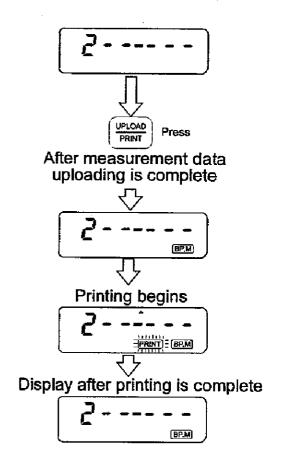

10.3.3. Time Setting

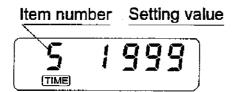

Explanation

- This setting sets the time in the blood pressure monitor (UA PC series).
- The display shown above appears when the TIME/ID key is pressed.
- The TIME indicator is displayed during time setting operation.
- An explanation of item numbers and setting values is shown in the table below.

Item number	Setting value range	Setting content
5	1997 ~ 2096	Year (4 digits)
6	1 ~ 12	Month
7	1 ~ 31	Day (setting range is automatically adjusted according to the year and month)
8	0 ~ 23	Hour (use 24-hour time)
9	00 ~ 59	Minutes (2 digits are always displayed)

Setting Procedure

- Step 1. When communication preparations are complete (refer to section 10.3.1 "Communication preparations"), press the TIME/ID key.
- Step 2. The TIME indicator and item number 5 will be displayed along with the setting value for the connected blood pressure monitor.
- Step 3. To change the setting value, press the
 ▲ ▼ keys.
 In this example, the setting value in item number 5 (year) will be changed to 1999.
- Step 4. When changing the item number, press the MENU key. When changing the setting value, press the ▲ ▼ keys.

 In this example, the setting value of item number 6 (month) will be changed to 4 (April).
- Step 5. When other changes (day, hour, minute) are necessary, select each item number and repeat the same basic procedure to change the setting.
- Step 6. When changes are complete press the ENTER key to confirm the new values.


 The new setting for the blood pressure monitor will be printed.

setting.


Please check the content of the

Step 7. If the CANCEL key is pressed without pressing the ENTER key, the time setting operation will be canceled and the setting of the blood pressure

monitor will not be changed.

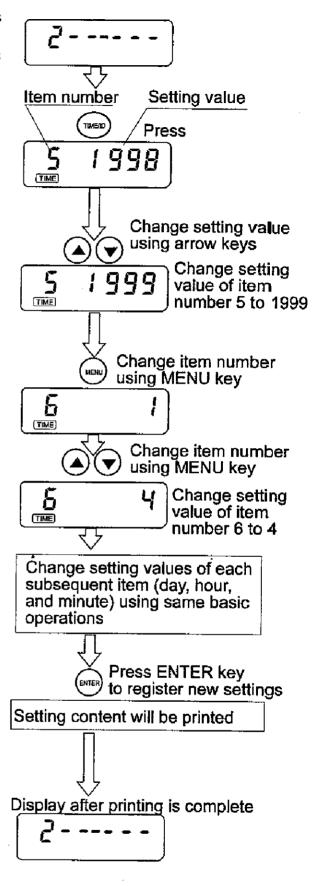

10.3.4. Alarm Setting for Measurement Time

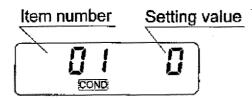

Explanation

- This function uses for the alarm at the measurement time. And it can set alarms up to six times.
- The above display appears when the INTERVAL key is pressed.
- The COND indicator is displayed during the operation.
- An explanation of the item numbers is provided in the table below.

Item number	Setting range	Setting conten	it
01	OFF, 0 ~ 23	First measurement time	(hour)
02	0 ~ 59	First measurement time	(minute)
03	OFF, 0 ~ 23	Second measurement time	(hour)
04	0 ~ 59	Second measurement time	(minute)
05	OFF, 0 ~ 23	Third measurement time	(hour)
06	0 ~ 59	Third measurement time	(minute)
07	OFF, 0 ~ 23	Fourth measurement time	(hour) ·
08	0 ~ 59	Fourth measurement time	(minute)
09	OFF, 0 ~ 23	Fifth measurement time	(hour)
10	0 ~ 59	Fifth measurement time	(minute)
11	OFF, 0 ~ 23	Sixth measurement time	(hour)
12	0 ~ 59	Sixth measurement time	(minute)


Refer to the operation manual of the blood pressure monitor (UA PC series) for details regarding this alarm function.

Setting Procedure

- Step 1. When communication preparations are complete (refer to section 10.3.1 "Communication preparations"), press the INTERVAL key.
- Step 2. The COND indicator will light and the time for the first alarm, set in the connected blood pressure monitor, will be displayed.
- Step 3. To change the setting value, press the ▲ ▼ keys.

 In this example the time for item number 1 (first alarm time (hour)) will be set to 10 (10:00AM).
- Step 4. Press the MENU key to change item numbers and press the ▲ ▼ keys to change setting values.
 In this example, the setting value in item number 2 (first alarm time (minutes)) will be set to 30 (30 minutes).
- Step 5. When changes in the second alarm and subsequent alarm are required, use the same procedure to select item numbers and set the hour and minute.
- Step 6. Press the ENTER key to register changes and exit the setting mode.

 The new setting for the blood pressure monitor will be printed. Please check the content of the setting.
- Step 7. If the CANCEL key is pressed without pressing the ENTER key, the current setting will be canceled. Then the blood pressure monitor will not be changed.

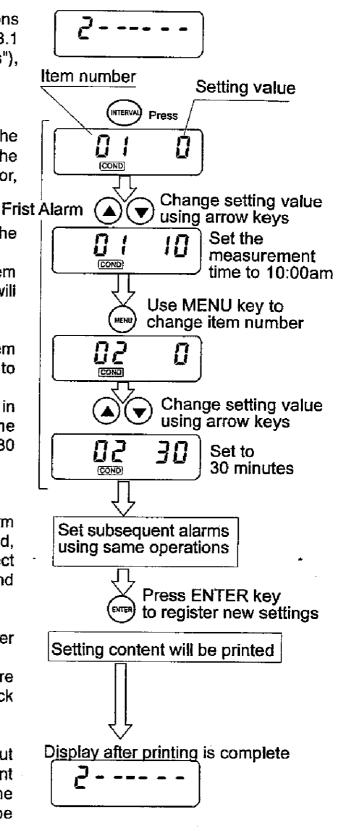

10.3.5. Memory Clear


Ex planation

The MEMORY CLEAR key is used to delete measurement data in the blood pressure monitor (UA PC series).

Setting Procedure

Step 1. When communication preparations are complete (refer to section 10.3.1 "Communication preparations"), press and hold the MEMORY CLEAR key. The buzzer will beep every second, and after 5 seconds a long buzzer sound will be heard. After measurement data has been cleared, the message "MEMORY CLEARED" will be printed. When measurement data has not been cleared. the message "MEMORY CLEAR FAILED" will be printed.

Step 2. When the MEMORY CLEAR key is not pressed and held for 5 seconds, the measurement data will not be deleted. In this case no message will be printed.

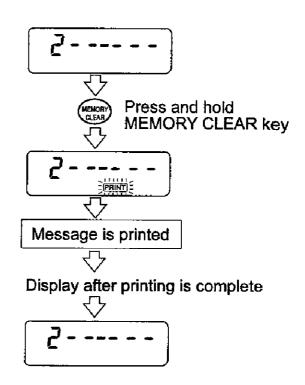

⚠ Caution!

Be careful when using these operations since measurement data cleared from the blood pressure monitor cannot be recovered. Be sure to execute measurement data processing before clearing the memory.

10.3.6. Printing of Setting List

Explanation


The alarm time set in the blood pressure monitor (UA PC series) will be printed.

Procedure


Step 1. When communication preparations are complete (refer to section 10.3.1 "Communication preparations"), press the ENTER key. Printing will begin.

Notice

When no blood pressure monitor is connected, no printing will be done.

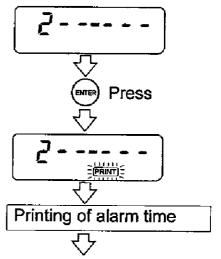


Display after printing is complete.

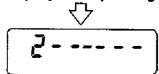

10.4. Model 3, Operation and Setting

When model 3 is selected (10.1.1 "Model selection") , operations are outlined in the following flowchart.

Model name: UB PC series

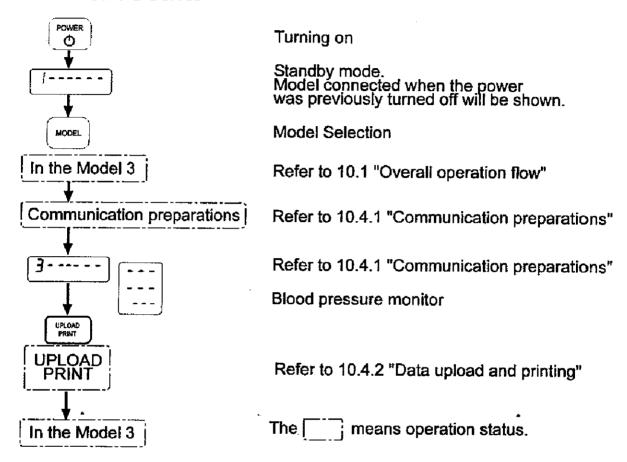

When model 3 is selected the functions that can be used with the printer terminal are limited. The functions and various setting that can be used are shown in the table below. For a detailed description of the setting methods, refer to the reference section listed in the table below.

Function name	Content	Reference section
UPLOAD/PRINT	Upload and print measurement data	10.4.2 "Data Upload Printing"

Only one-way communication functions can be used with model 3. The blood pressure monitor time setting and measurement data clearing operation are done on the blood pressure monitor side. Use the printer after reading the operation manual of the blood pressure monitor (UB PC series).


A special communications adapter, sold separately, is required with model 3.

⚠ Caution!

When the unit is left with the communication cable inserted in the blood pressure monitor, the batteries will rapidly drain. Therefore, when communications are complete, unplug the cable. When the cable is unplugged the blood pressure monitor will return to the clock display.

10.4.1. Communication Preparations

- Step 1. Connect the special communications adapter, sold separately, to the connector on the back of the printer terminal (refer to section 9. "Preparation").
- Step 2. Connect the mini-pin jack end of the special communications adapter to the blood pressure monitor (UB PC series).
- Step 3. The blood pressure monitor display will change to [---
- Step 4. When the display changes communication preparations are complete.

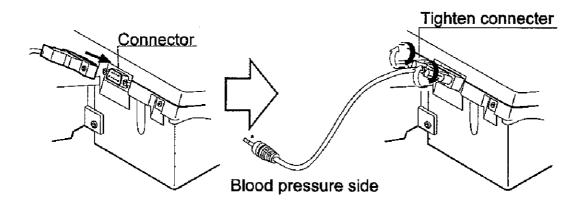

⚠ Caution!

- When the unit is left with the communication cable inserted in the blood pressure monitor, the batteries will rapidly drain. Therefore, when communications are complete, remove the cable. When the cable is removed the blood pressure monitor will return to the clock display.
- The preparation procedure is different about several blood pressure monitors. Refer to each instruction manual of the blood pressure monitor.

10.4.2. Data Upload and Printing

Explanation

 Press the UPLOAD PRINT key then immediately press the start switch on the blood pressure monitor. Measurement data will be uploaded from the blood pressure monitor. Measurement results will then be printed using the standard report format.


With model 3 the report format cannot be selected.

- When measurement data is uploaded, the BPM indicator is displayed Measurement data uploading is complete when this indicator is displayed. The communication cable can be unplugged from the blood pressure monitor at this time.

Procedure

- Step 1. Confirm that communication preparations are complete (refer to section 10.4.1 "Communication preparations").
- Step 2. Press the UPLOAD/PRINT key. The printer will enter the communication wait mode.
- Step 3. Immediately press the start switch of the blood pressure monitor.

 Measurement data will be sent from the blood pressure monitor.
- Step 4. The measurement data will be uploaded to the printer and the BP.M indicator is displayed when uploading is complete.
- Step 5. Next, printing will begin. The PRINT indicator will blink during printing.
- Step 6. The PRINT indicator will go out when printing is complete.

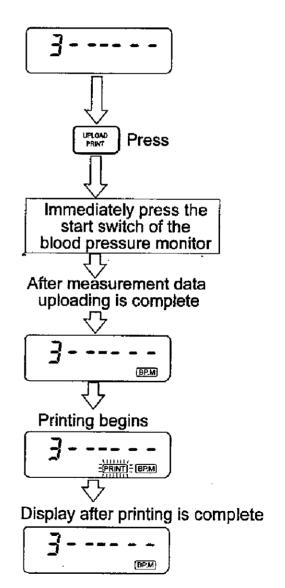

11. Maintenance

11.1. Cleaning

- When cleaning the unit, first turn the power off and unplug the AC adapter from the printer terminal. In addition, remove the batteries, if present, after the power is turned off.
- Wipe the unit with a dry cloth.
- When the printer is too dirty, clean using a well-wrung cloth wet with a neutral cleaner.
- The printer terminal is not waterproof. Do not wet the unit when cleaning.
- Do not clean using thinner or other organic solvents.

11.2. Before Requesting Repairs

- Do not open the case and do not allow anyone other than authorized personnel to repair the unit. Such repairs may cause damage or fires.
- Before assuming that the unit is broken, try solving the problem by matching the
 error number displayed on the LCD with the measures described below. If the
 problem persists, contact the nearest business representative of A & D Co., Ltd.

Notice

When an error is displayed the power must first be turned off before operation can be started again. Turn the power off, then on again.

Error number	Error content	Cause or countermeasure for problem
E01		 Remove the AC adapter or batteries and insert again. Turn the power on.
E02	Memory error	If the error re-appears, send the unit in for repair.
E11	Communication error	 Is the communication cable properly connected? Do the connected blood pressure monitor and the selected model match? Is the communication method correct? (when using model 1) Have communication preparations been properly carried out? (preparations are different depending on the model)
E12	Infrared communication device error	 Remove the AC adapter or batteries and insert again. Turn the power on. If the error re-appears, send the unit in for repair.
E21	Paper jam	- Is the paper jammed?


12. Options


Part name	Order code
Printer paper (5 rolls)	AX-PP147-S
Communication cable (for models 1 and 2) The same cable as the accessory cable	AX-KO1502


13. Appendix

13.1. External Diagrams


TM-2480

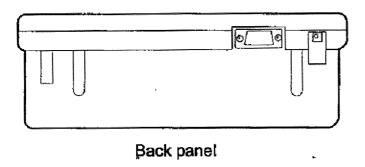


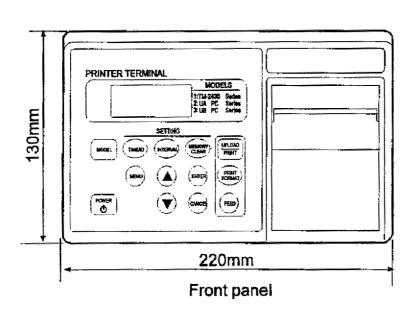


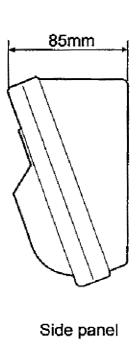

TM-2481

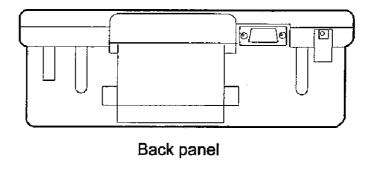


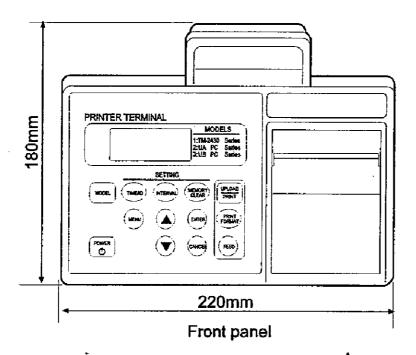


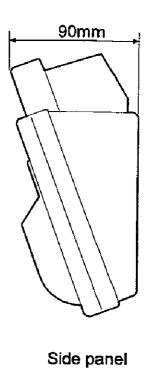



13.2. Printout Examples

Printout example for TM-2430 series (printing is done in the sequence shown below every 24 hours)

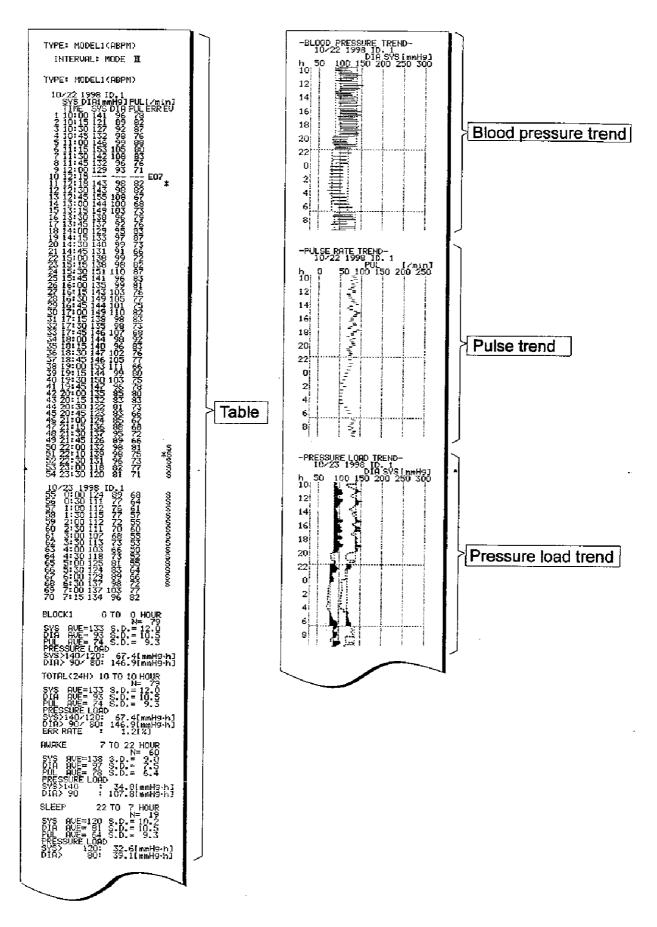

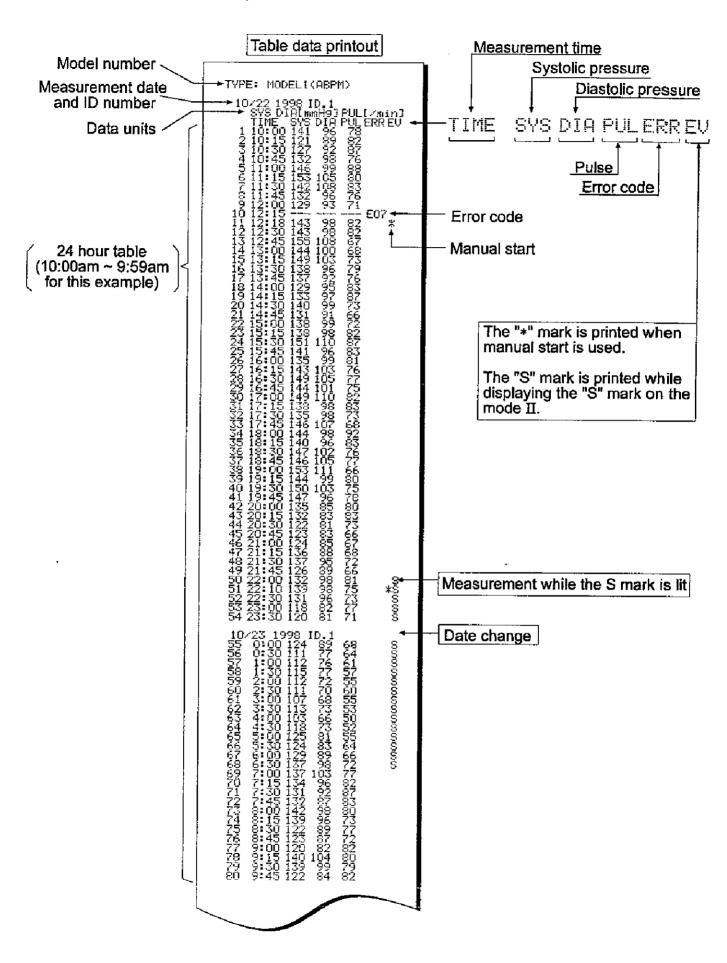

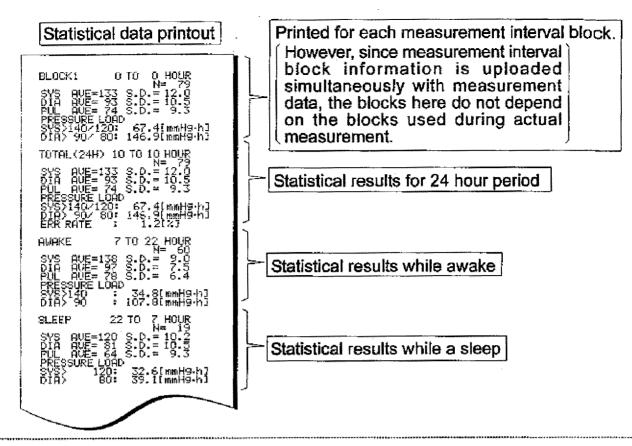

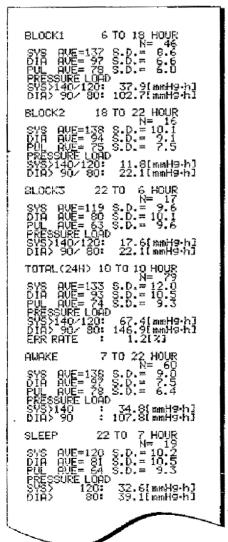

This example shows the printout with only the measurement interval blocks changed. The base data used to obtain the above statistical data was also used here.


Changed to 3 blocks


The method for changing the measurement interval into 3 blocks, as shown in the example at right, is shown below.


Change method

- Step 1. Set the blood pressure monitor to the connected mode.
- Step 2. Change the measurement interval blocks of the blood pressure monitor (refer to section 10.2.5 "Measurement interval settings")
- Step 3. Press the UPLOAD PRINT key.


The TOTAL (24H)/AWAKE/SLEEP content will be the same as the above printout.

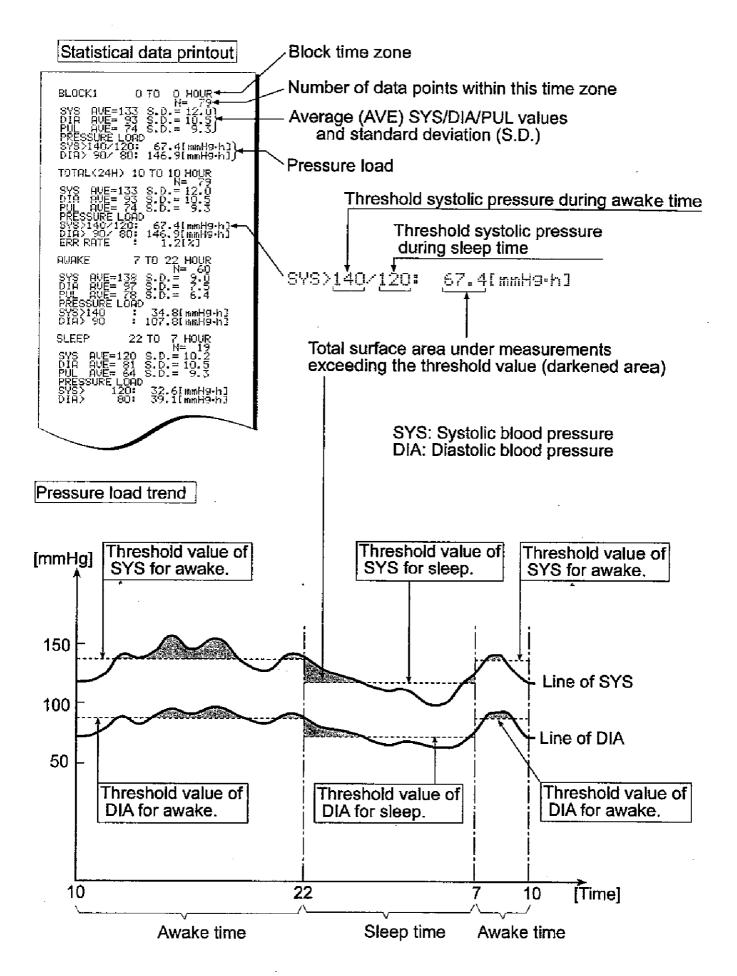


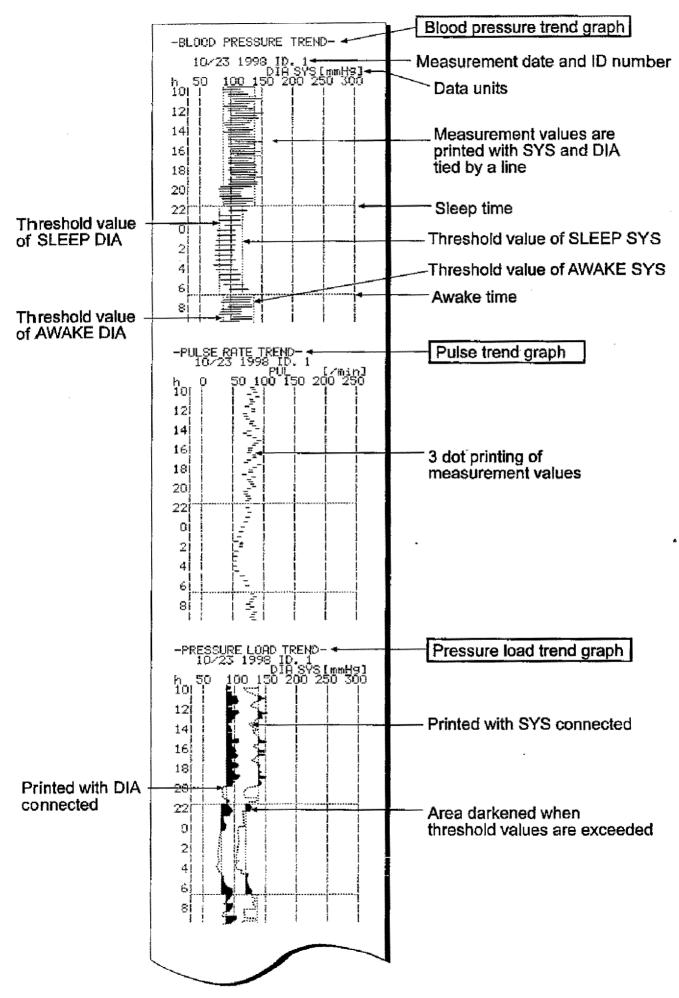


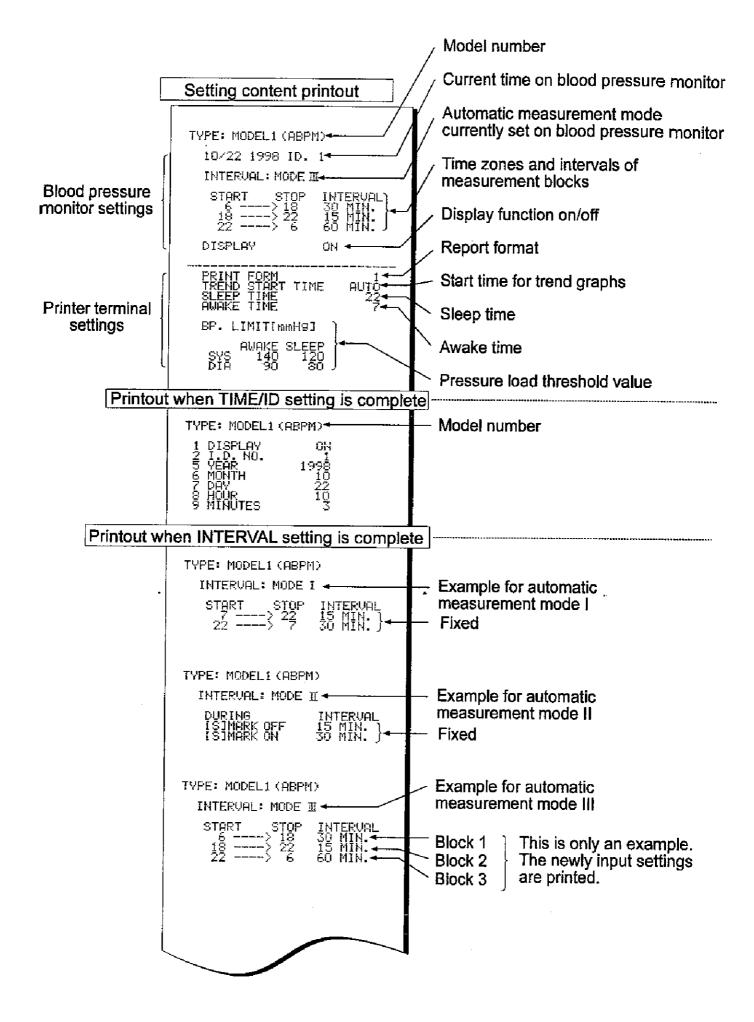


Printout example for UA PC series /UB PC series

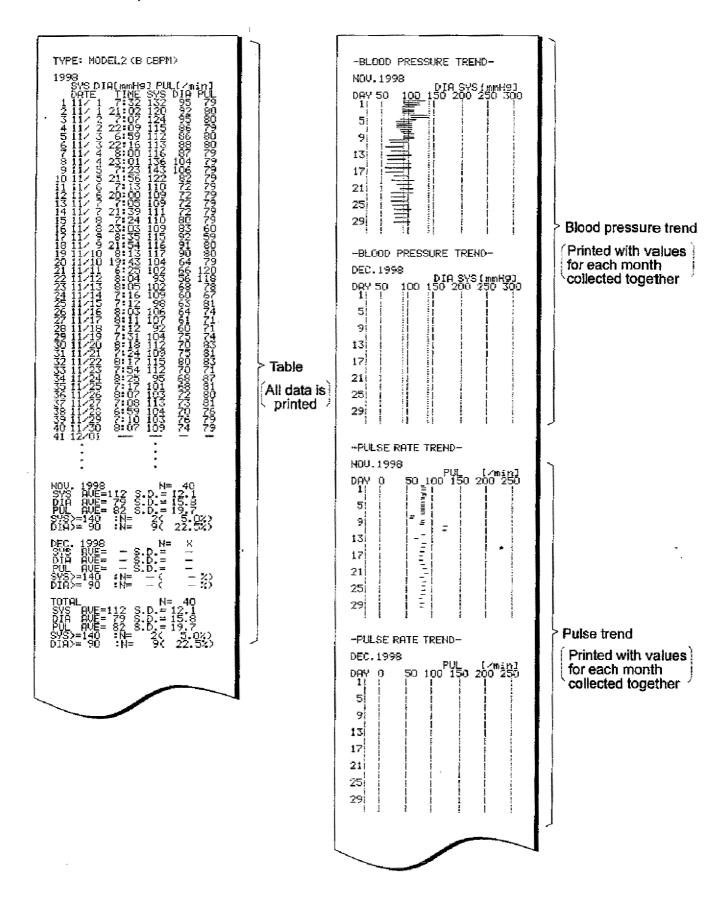

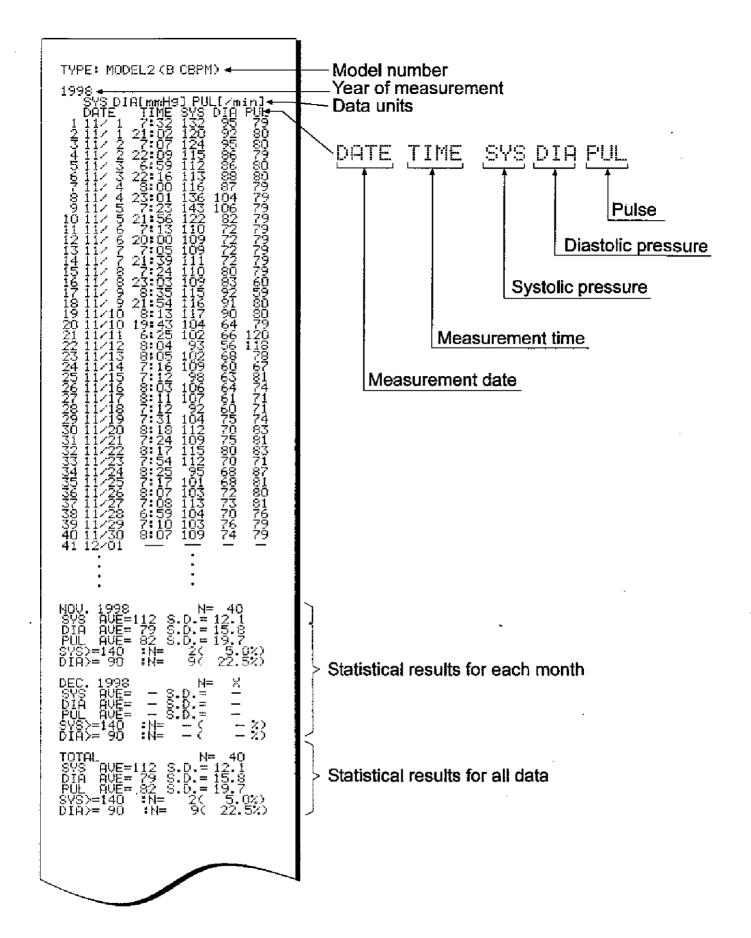

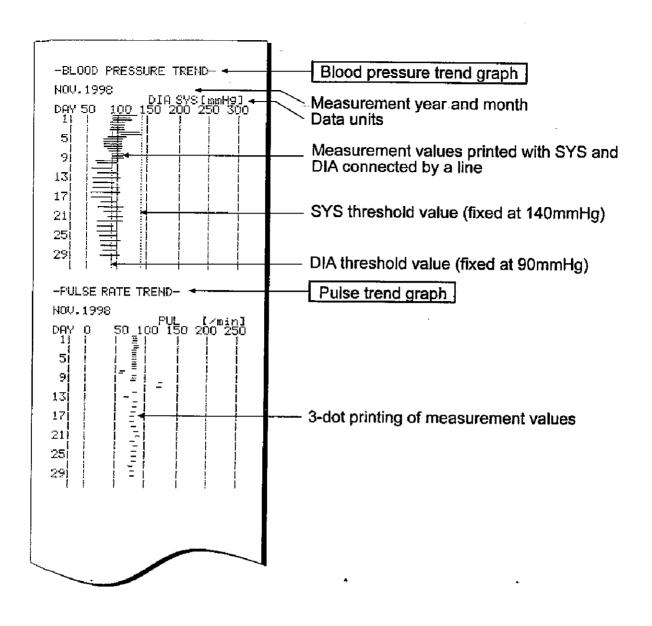

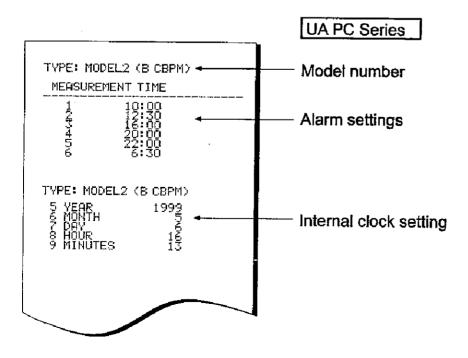



A&D Co., LTD. 1-243 Asahi, Kitamoto-shi, Saitama, 364-8585 JAPAN Telephone: [81] (485) 93-1111 Fax: [81] (485) 93-1119

A&D INSTRUMENTS LTD.

<Authorized Representative Established in the European Community>Abingdon Science Park, Abingdon, Oxford OX14 3YS ENGLAND Telephone: [44] (1235) 550420 Fax: [44] (1235) 550485