

DIGITAL BLOOD PRESSURE MONITOR

UA-787Plus

Clinically
Validated*

IHB



*Validated as UA-787

The Monitor

- Smooth and Fast Measurement
- Irregular Heart Beat (IHB) Indicator
- Blood Pressure Classification Indicator
- Large 60-Memory + Average Reading including Time and Date
- Newly designed e-Cuff for Easy Fastening and Comfortable Measurement
- Large Comprehensive 3-line Display
- Clock Display and Reminder Alarm



CERTIFIED



CE 0366

A&D
A&D Medical

Technology You Can Trust

<http://www.aandj.jp>

DIGITAL BLOOD PRESSURE MONITOR UA-787Plus

Smooth and Fast Measurement

Patients may feel a little uncomfortable during blood pressure measurement because of the pressure applied to the arm and the noise while inflating the cuff. Our newest technology has made those feelings of discomfort a thing of the past. The UA-787 Plus requires a shorter measurement time and its newly developed semi-hard shell cuff also caters for a wider range of arm sizes. The new whisper silent pump further aids smooth and fast measurement.

Irregular Heart Beat (IHB) Indicator

This lights up if an IHB is detected during a measurement. The indicator alerts users of an irregular heartbeat, which could be linked with cardiovascular disease. (Note: This is for warning purposes only. A&D recommends contacting a physician if the user observes this indicator frequently.)

Blood Pressure Classification Indicator

The World Health Organization (WHO) distinguishes 6 different blood pressure stages. This monitor automatically indicates your classification (Optimal to Severe hypertension) after each measurement is completed.

Large 60-Memory + Average Reading

The UA-787 Plus has a 60 reading memory capacity including time & date for monitoring the trend of your measurements. An average of all the stored data is shown at each start up, offering the user their comprehensive measurement history at a glance.

Newly designed e-Cuff

A newly designed semi-hard shell cuff promises easy fastening around your arm as well as comfortable blood pressure measurements for patients in a wider arm range (17 to 32cm).

Large Comprehensive 3-line Display


The UA-787 Plus is equipped with a comprehensive LCD display so you can feel confident about your blood pressure measurement and its results. A detailed result readout (Systolic, Diastolic & pulse) is shown simultaneously. Additional information such as the inflation level, WHO classification & IHB, as well as the SYS/DIA readout, offer the user a wide range of details.

Clock Display with Reminder Alarm

The Current Time is displayed when the unit is not in use. You can set the alarm for 3 separate times during the day to remind yourself of your daily measurements.

WHO blood pressure classification

Standards for assessment of high blood pressure, without regard to age, have been established by the World Health Organization (WHO), as shown in the chart below.

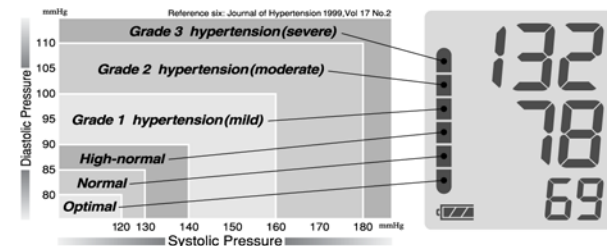

A&D
A&D Medical
Tokyo JAPAN
<http://www.aandj.jp>

Technology You Can Trust

*UA787PlusCA-ADCC-01-PR1-05a10
REF*CARD-UA787Plus